
phoenix

implante dental de conexión externa
external connection dental implant

www.ilerimplant.com

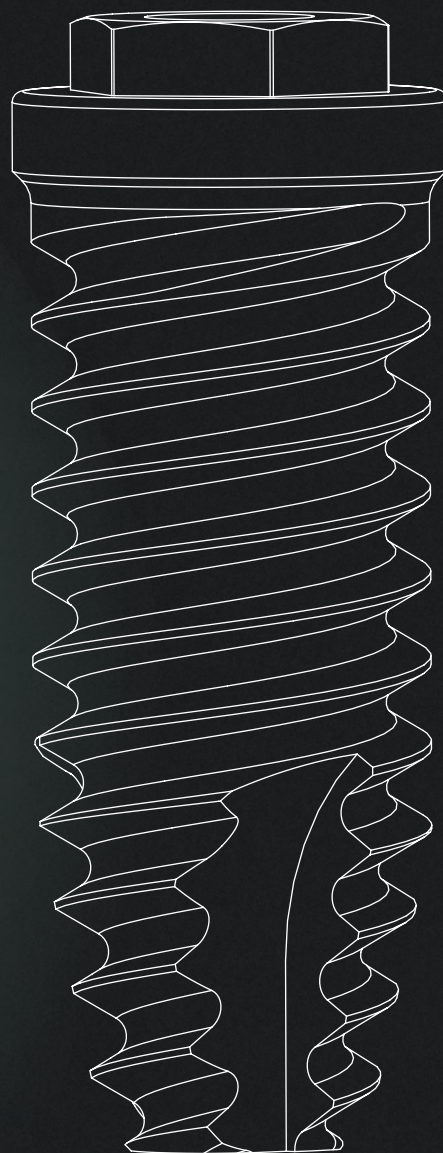
Made in Spain

phoenix

implante dental de conexión externa
external connection dental implant

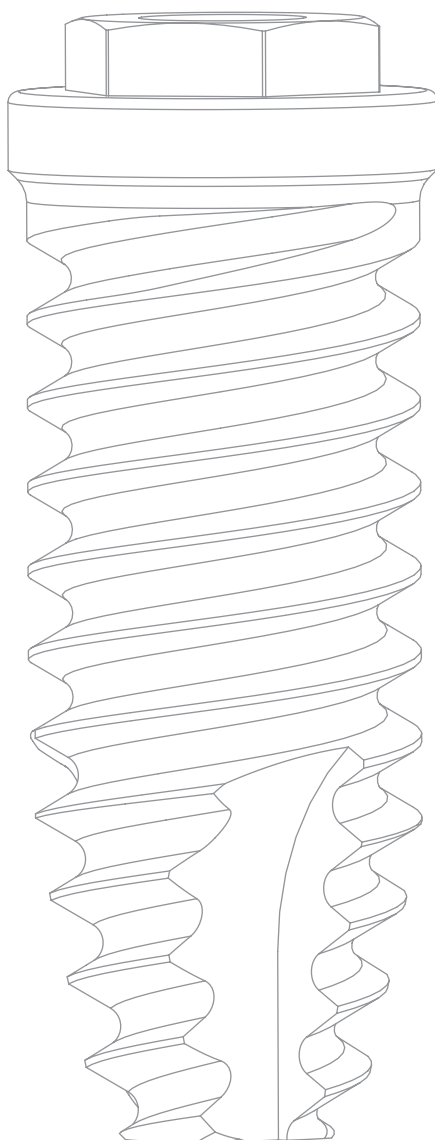Sistema de
implantes
Implant system

Índice - Summary

SISTEMA DE IMPLANTES

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ATTACHMENTS

- gmi phoenix healing abutments 14
- gmi phoenix esthetic angled abutments 16
- gmi phoenix straight abutments 17
- gmi phoenix millable-implant carrier abutments 18
- gmi phoenix ball abutments 19
- gmi phoenix multi-esthetic abutments 21
- gmi phoenix replicas 24
- gmi phoenix impression coping 25
- gmi phoenix castables 26
- gmi phoenix clinic and laboratory screws 28
- Recommended tightening torques 29
- gmi phoenix surgical kit 30

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
gmi phoenix external connection implants are machined in CP grade IV titanium of the highest quality, thus ensuring a high chemical stability and excellent biocompatibility. The titanium oxide layer that forms on the implant surface after the passivation process, and the careful design of the outer surface of the thread, guarantee good osseointegration and high stability.

The external dual inlet thread of the implant helps to reduce surgery time. Moreover, the rounded shapes, self-threading millings and slight tapering facilitate implant insertion, reducing tension at the bone-implant interface and thus preventing threading stress-related problems.

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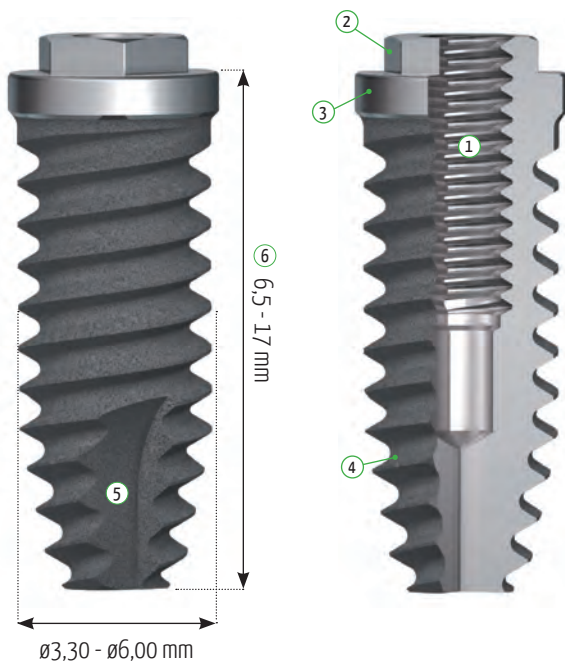- ② **External hexagon**
 Narrow: B/F 2.40 mm x 1 mm
 Standard: B/F 2.70 mm x 0.7 mm
 Wide: B/F 2.70 mm x 0.7 mm
 B/F 3.40 mm x 0.7 mm
 Extra wide: B/F 3.40 x 0.9 mm

- ③ **Platforms**
 1.60 mm pitch
 ADS surface treatment

- ⑤ **Self threading millings**
 Variable length
 3 millings at 120°

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


To promote the adhesion and growth of bone cells, the outer surface of the implant has been treated with the exclusive ADS® treatment (Advanced Double-Grip Surface) which combines a white corundum micro-bubble treatment and acid etching, for a non-smooth roughness that maximizes the contact surface between implant and bone and therefore improves stability and adhesion between them.

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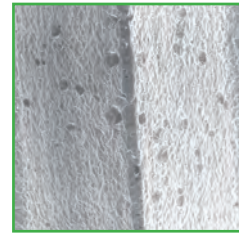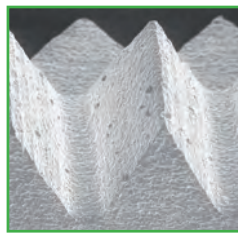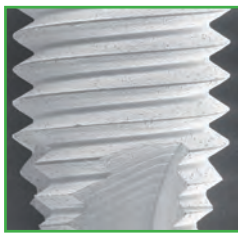
The ADS® treatment is the result of a series of collaborations with institutes and universities, as well as our own research and tests, that allow us to monitor and validate the optimal biological response to the implant surface: chemical (X-ray scattering studies to find out the composition of the implant and its surface) and topography studies (roughness studies) have been carried out, as well as biological (in vivo response studies in animals) and clinical studies.

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X-ray scattering studies show that the surface of the implants is only made up of titanium and oxygen, which are the basic components of the oxide layer that forms in the implant passivation process and that gives it its excellent corrosion resistance.

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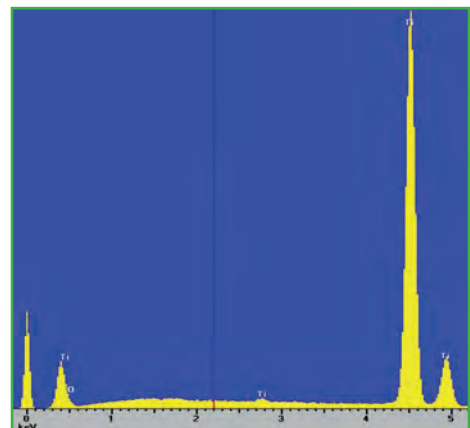
In order to control and identify the topography of the implant surface, roughness studies were performed with profilometer, as well as with techniques such as SEM (Scanning Electron Microscope) and CLSM (Confocal Laser Scanning Microscope), allowing to obtain and control the mean roughness (Ra) values within the parameters recommended in international publications.



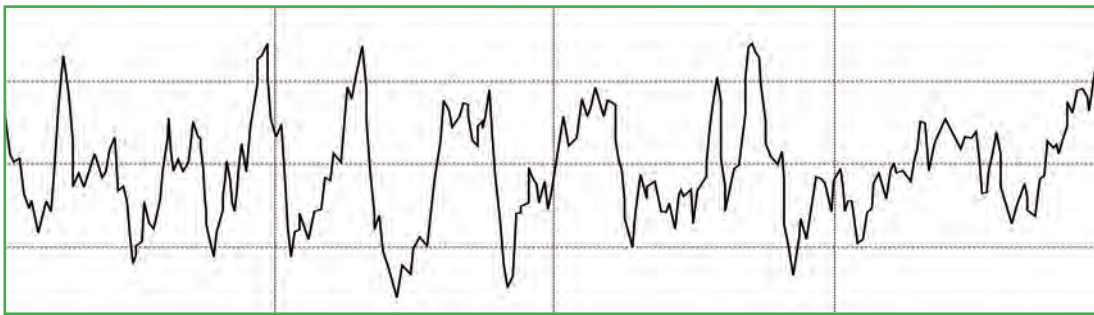
Perfil de rugosidad obtenido con HOMMEL-WERKE T1000 • Roughness profile obtained with HOMMEL-WERKE T1000.

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To determine the biological response of GMI implants, in vivo studies were conducted by inserting the implants on animals, leaving them without load during the healing process and performing a histological study to prove the excellent response of the bone cells and the complete osseointegration of the implant.

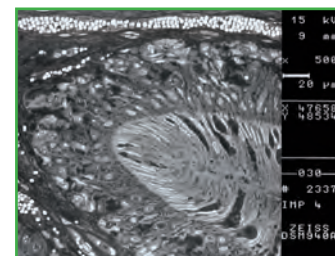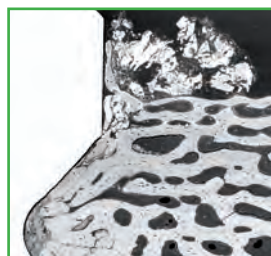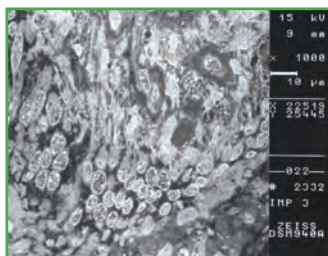SEM-BSE: detalles ultraestructurales de las células en biofilm sobre el implante. • SEM-BSE: ultrastructural details of cells in biofilm on the implant.

Protocolo de inserción de implantes gmi phoenix - Insertion sequence for gmi phoenix implants

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gmi phoenix external connection implants come in a screen-printed cardboard box, with the instructions for use, traceability labels and a double packaging containing the implant system, with the primary and secondary Pyrex glass packaging to ensure the inert performance of the container. Once packaged in a clean room, they are sterilized with gamma radiation according to regulations.

The new packaging system of the **gmi phoenix** implants consists of a titanium support which only houses the implant-implant carrier set, prepared to be inserted directly into the mouth, using a single hex-3.00 mm implant carrier wrench, designed for this purpose. This new system offers the odontologist the following advantages: more control and a better view of the implant insertion process, it is easier to use in reduced interdental spaces and a reduction in surgery time.

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1. Check sticker indicating implant diameter and length and open the tab on the carton.
2. Remove the secondary packaging, instructions for use and product identification labels.
3. Check the integrity of the safety seal and unscrew the plug on the secondary packaging.
4. Remove the primary packaging from the secondary packaging avoiding hitting it against a hard surface.
5. Extract and save the primary packaging plug making a lateral movement.
6. Insert the HEX-3.00 mm implant carrier wrench until a slight retention is noticed, addressing the notches on the key with the faces of the implant hex broaching.
7. Check that the key is fully inserted and turn slightly while gently pulling on the assembly.
8. Place the implant in the bone bed by applying a maximum torque of 35 N-cm until the treated area is at the crestal or subcrestal level. Remove the key.
9. Unscrew the clinic screw and dismount the implant carrier by turning slightly if necessary.
10. Remove the cover screw from the primary packaging plug, using the HEX-0.80 mm ratchet wrench for Ø3.30 platform and HEX-1.20 mm for the other platforms.
11. Screw the cover screw to the implant by applying a maximum torque of 15 N-cm. Remove key and suture incision.

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

Following is a description of the drilling sequences for the different external connection implant models, as well as the recommended drilling conditions for their use.

- Lance-shaped: 1200 - 1500 rpm
- Pilot drill: 900 - 1200 rpm
- Guide drill: 800 rpm
- Twist drills:
 - Ø2.80 mm → 500 - 700 rpm.
 - Ø3.00 and Ø3.50 mm → 400 - 700 rpm.
 - Ø4.25 mm → 400 - 600 rpm.
 - Ø5.10 and Ø5.40 mm → 300 - 500 rpm.
- For drilling, external cooling with saline solution is required.

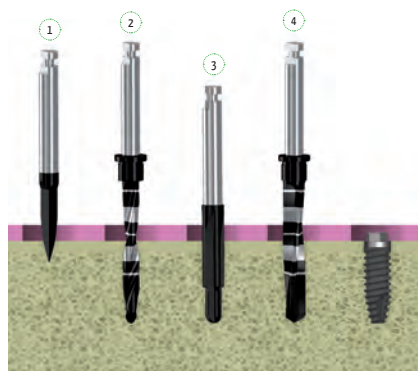
phoenix Ø3,30		
	Descr.	Ref.
1 Fresa lanceolada	Lance-shaped drill	KYF OC 1224
2 Fresa piloto	Pilot drill	KYF OC 2221
3 Fresa guía	Guide drill	KYF OC 3227
4 Fresa Ø2,80 mm	Twist drill Ø2.80 mm	KYF OC 0128

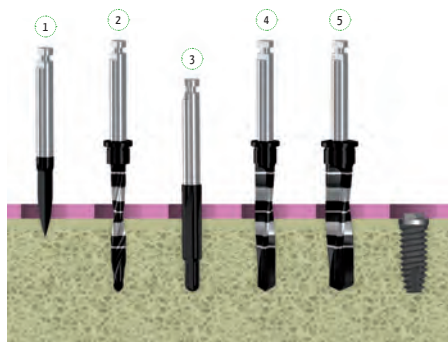
phoenix Ø3,75		
	Descr.	Ref.
1 Fresa lanceolada	Lance-shaped drill	KYF OC 1224
2 Fresa piloto	Pilot drill	KYF OC 2221
3 Fresa guía	Guide drill	KYF OC 3227
4 Fresa Ø2,80 mm	Twist drill Ø2.80 mm	KYF OC 0128
5 Fresa Ø3,25 mm	Twist drill Ø3.25 mm	KYF OC 0132

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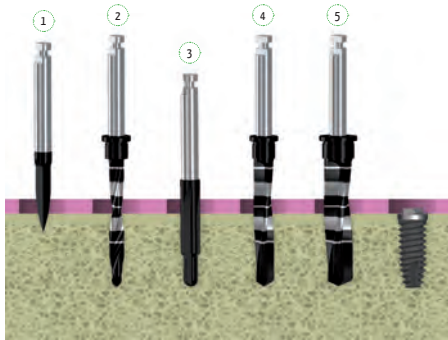
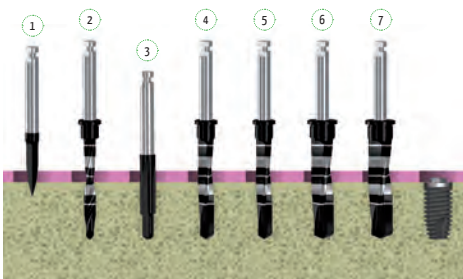
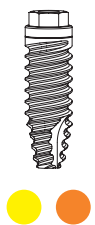
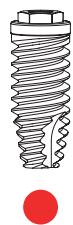
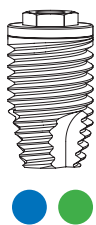
The **gmi phoenix** implant range was designed and tested to be inserted as a single restoration, with a 30-degree maximum inclination to the perpendicular plane to the occlusal plane, in accordance with the following diagram.

	 Ø3,30-Ø3,75 mm Incisivos (Excepto centrales maxilares) Incisors (Except central maxillary incisors)	 Ø4,00 mm Caninos • Canines Premolares • Premolars Incisivos centrales maxilares Central maxillary incisors	 Ø5,00-Ø6,00 mm Molares Molars
---	---	--	---

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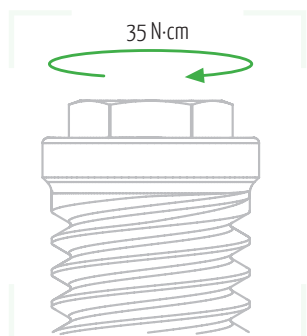
The GMI PHOENIX implant range is designed to be inserted applying a 35-40 N-cm torque, although it can bear a maximum tightening torque of 60 N-cm. GMI recommends manual insertion of the implant with a dynamometric ratchet wrench to ensure that these values are not exceeded.

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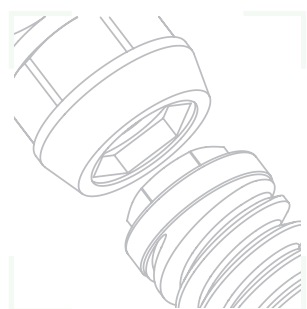The GMI PHOENIX range of implants and prosthetic attachments are designed and manufactured with tolerances intended to maintain adjustments that minimize the clearance between components and therefore the stability of the connection. GMI recommends the use of original attachments to ensure a perfect fit between the components of the restoration.

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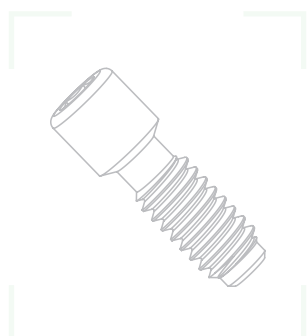
To prevent prior deformation of the clinic screw, GMI recommends restricting its use only for the final tightening of the prosthesis. To make adjustments during the prosthetic procedures, use of the laboratory screw is recommended.

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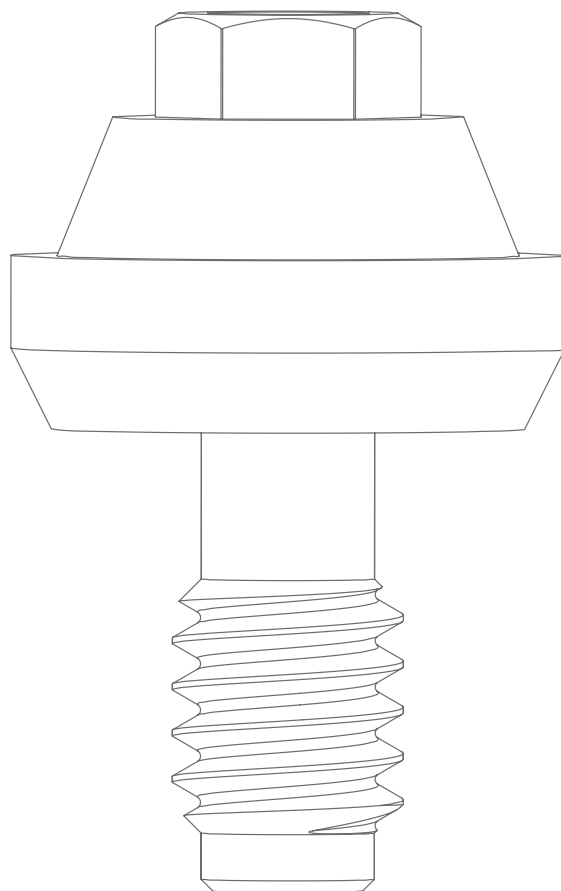Once the repair phase of the tissue supporting the implant has finished, in the case of implants fully or partially covered by soft tissue there must be a mucosal channel or tunnel enabling connection of the implant to the prosthesis or secondary structure. Healing abutments are used to form the mucosal route, once the soft tissue has been opened by incision and the cover screw has been withdrawn, by screwing them to the implant.

Key features

- Grade V titanium machining as per ISO-5832 standard.
- Tightening with b/f 1.20 mm hex wrench.
- Alternative tightening with a flat tip screwdriver.
- Available in emergency heights from 2.00 to 5.00 mm.
- Slightly tapered in the head to facilitate extraction.
- Apical bit to facilitate the start of threading.

Pilares de cicatrización para plataforma Ø3,30 mm Healing abutments for Ø3.30 mm platform

Hex. impl.	Altura (mm) Height (mm)	Ref.
◇ 2,40 mm	2,00	KDP OC 0023
	3,00	KDP OC 0024
	4,00	KDP OC 0025
	5,00	KDP OC 0026

Pilares de cicatrización para plataforma Ø4,10 mm Healing abutments for Ø4.10 mm platform

Hex. impl.	Altura (mm) Height (mm)	Ref.
◇ 2,70 mm	2,00	KDP OC 0005
	3,00	KDP OC 0006
	4,00	KDP OC 0007
	5,00	KDP OC 0008

Pilares de cicatrización para plataforma Ø5,10 mm
Healing abutments for Ø5.10 mm platform

Hex. impl.	Altura (mm) Height (mm)	Ref.
○ 2,70 mm	2,00	KDP OC 0011
	3,00	KDP OC 0013
	4,00	KDP OC 0015
	5,00	KDP OC 0017
○ 3,40 mm	2,00	KDP OC 0012
	3,00	KDP OC 0014
	4,00	KDP OC 0016
	5,00	KDP OC 0018

Pilares de cicatrización para plataforma Ø6,00 mm
Healing abutments for Ø6.00 mm platform

Hex. impl.	Altura (mm) Height (mm)	Ref.
○ 3,40 mm	2,00	KDP OC 0029
	3,00	KDP OC 0030
	4,00	KDP OC 0031
	5,00	KDP OC 0032

Pilares inclinados estéticos gmi phoenix - gmi phoenix esthetic angled abutments

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When it is necessary to correct an extreme angular position of the implant inside the bone structure in relation to adjacent natural or prosthetic teeth, angled abutments should be used, that once mounted and oriented on the implant will allow the parallelization of the new crown with other structures.

Key features

- Grade V titanium machining as per ISO-5832 standard.
- Allow non-parallelism correction up to 20 degrees.
- 12-position interior for increased flexibility.
- Clinic screw included.
- Abutment shoulder adapted to the shape of the gum.

Pilares inclinados estéticos para plataforma Ø3,30 mm Esthetic angled abutments for Ø3.30 mm platform

Hex. impl.	Tipo conexión Connection type	Ang. (°)	Ref.
⬡ 2,40 mm	Antir. • Non-rot.	15°	KDM of 0060
	Antir. • Non-rot.	20°	KDM of 0069

Pilares inclinados estéticos para plataforma Ø4,10 mm Esthetic angled abutments for Ø4.10 mm platform

Hex. impl.	Tipo conexión Connection type	Ang. (°)	Ref.
⬡ 2,70 mm	Antir. • Non-rot.	15°	KDM of 0002
	Antir. • Non-rot.	20°	KDM of 0070

Pilares inclinados estéticos para plataforma Ø5,10 mm Esthetic angled abutments for Ø5.10 mm platform

Hex. impl.	Tipo conexión Connection type	Ang. (°)	Ref.
⬡ 2,70 mm	Antir. • Non-rot.	15°	KDM of 0008
	Antir. • Non-rot.	20°	KDM of 0071
⬡ 3,40 mm	Antir. • Non-rot.	15°	KDM of 0014
	Antir. • Non-rot.	20°	KDM of 0072

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Once the mucosal tunnel has been formed by the healing abutment and after removing the cap from this location, it is appropriate to use the abutment, which is defined as the machined tapered attachment that, after being fixed with the clinic screw to the top of the implant, acts as a support for the clinic crown and is attached to it using cement.

Key features

- Grade V titanium machining as per ISO-5832 standard.
- Anti-rotational planes to increase retention.
- Straight abutments in hex rotary and non-rotary version.
- Clinic screw included.

Pilares rectos para plataforma ø3,30 mm Straight abutments for ø3.30 mm platform

Hex. impl.	Tipo conexión Connection type	Ref.
⬡ 2,40 mm	Antir. • Non-rot.	KDM OF 0058

Pilares rectos para plataforma ø4,10 mm Straight abutments for ø4.10 mm platform

Hex. impl.	Tipo conexión Connection type	Ref.
⬡ 2,70 mm	Antir. • Non-rot.	KDM OF 0034

Pilares rectos para plataforma ø5,10 mm Straight abutments for ø5.10 mm platform

Hex. impl.	Tipo conexión Connection type	Ref.
⬡ 2,70 mm	Antir. • Non-rot.	KDM OF 0035
⬡ 3,40 mm	Antir. • Non-rot.	KDM OF 0036

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When the use of a custom abutment is required, millable abutments can be used. By cutting them in any direction, they can be suited to the needs of the prosthetic. In addition and as a complementary function, millable transfer abutments can be used like impression coping, threading them to the implant screw through the short impression screw and using them to make the negative of the patient's mouth for subsequent positivation on plaster.

Key features

- Grade V titanium machining as per ISO-5832 standard.
- They allow cutting in any direction to suit the needs of the prosthetic.
- Compatible with b/f 4.00 mm hex wrench for use as a implant carrier.
- Available only in non-rotary version.
- Possibility to use as impression coping with short impression screw.
- Clinic screw included.

Pilares fresables para plataformas ø3,30 / 4,10 / 5,10 / 6,00 mm Millable abutments for ø3.30 / 4.10 / 5.10 / 6.00 mm platforms

Hex. impl.	Plataforma Platform	Ref.
⊕ 2,40 mm	ø3.30	KDB OF 0004
⊕ 2,70 mm	ø4.10	KDB OF 0001
	ø5.10	KDB OF 0002
	ø5.10	KDB OF 0003
⊕ 3,40 mm	ø6.00	KDB OF 0005

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Once the healing abutment has formed the mucosal tunnel and the healing abutment has been removed from this location, the corresponding ball abutment is used, which is defined as a machined attachment being fixed to the top of the implant and serving as the retention element for the removable prosthesis, with the previous addition of the retention caps that will be fixed to the top spherical part of the ball abutment.

Key features

- Grade V titanium machining as per ISO-5832 standard.
- Available in emergency heights from 1.50 to 3.00 mm.
- Pink retention cap and metal container included.
- Protection disk included.
- Tightening with b/f 2.45 mm hex wrench, part no.: KYL OC 0065.
- Instruments for caps: Universal handle (RH124ICP), Insertion tool for caps (RH185IAC) and Removal tool for caps (RH091EC).

Pilares bola para plataforma ø4,10 mm Ball abutment for ø4.10 mm platform		
Hex. impl.	Altura (mm) Height (mm)	Ref.
○ 2,70 mm	3,00	KDF OF 0045
	2,50	KDF OF 0046
	2,00	KDF OF 0047
	1,50	KDF OF 0048

Pilares bola para plataforma ø5,10 mm Ball abutment for ø5.10 mm platform		
Hex. impl.	Altura (mm) Height (mm)	Ref.
○ 3,40 mm	3,00	KDF OF 0049
	2,50	KDF OF 0050
	2,00	KDF OF 0051
	1,50	KDF OF 0052
○ 2,70 mm	3,00	KDF OF 0053
	2,50	KDF OF 0054
	2,00	KDF OF 0055
	1,50	KDF OF 0056

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Following the first stage of surgery and therefore after implant osseointegration is achieved, the multi-esthetic abutment (MEA) is used, that is defined as the machined attachment that, after being attached to the crown of the implant, serves as a support for a multiple restoration, thus decreasing the height between the implant connection and the structure connection, and providing a passive adjustment of the structure.

Key features

- Grade V titanium machining as per ISO-5832 standard.
- Available in emergency heights from 1.50 to 4.50 mm.
- One-piece straight abutments in rotary model.
- Tightened with hex keys:
PL-Ø3.30 and Ø4.10: B/F 2.00 mm (KYL OC 0097)
PL-Ø5.10: B/F 2.50 mm (KYL OF 0098)
- Impression coping with screw included.
- Castable with clinic screw included.
- Temporary cylinder in grade V titanium with clinic screw included.
- Angled multi-esthetic abutments allow non-parallelism up to 30°.

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Once the transfer of the position has been made using the impression coping from the biological medium to the laboratory working model, it is necessary to use replicas, which are defined as the machined attachments designed on one hand to reproduce the implant position in the mouth of the patient to a working model, and on the other to serve as a connection model for the construction in the laboratory of the structure to carry out the restoration.

Key features

- Machined replicas in high quality stainless steel.
- Anti-rotation millings and vertical retention neck for a good fixation.

Réplicas de implante para plataformas ø3,30 / 4,10 / 5,10 / 6,00 mm Implant replicas for ø3.30 / 4.10 / 5.10 / 6.00 mm platforms

Hex. impl.	Plataforma (mm) Platform (mm)	Ref.
○ 2,40 mm	ø3,30	KDR oC 0010
○ 2,70 mm	ø4,10	KDR oC 0002
	ø5,10	KDR oC 0006
○ 3,40 mm	ø5,10	KDR oC 0008
	ø6,00	KDR oC 0011

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Impression coping are defined as the machined attachments that, connecting to the coronal part of the implant within the oral cavity and axially fixed by the impression screw, is used to make the transfer of the position of the implant in the biological medium to a laboratory model, on which the prosthetic will work to create the dental crowns or prostheses. This process is achieved using impression materials placed in an appropriate manner, that will harden in the oral cavity and produce the negative of the patient's mouth. Once extracted, replicas are placed and the casting plaster is poured to obtain the positive mould where the replica will be placed in the original position of the implant in the mouth.

Key features

- Machined open tray impression coping in grade V titanium according to ISO-5832.
- Impression screws in grade V titanium according to ISO-5832.
- Tightening of impression screws with b/f 1.20 mm hex wrench.
- Apical bit to facilitate the start of threading.
- Anti-rotation millings and vertical retention neck for a good fixation.
- Impression coping screw included.
- Long impression screw with longitudinal millings for manual tightening.
- Closed tray impression coping machined in technical plastic.

Tomas impresión cubeta abierta (CA) plataforma Ø3,30 mm Open tray (OT) impression coping for Ø3.30 mm platform

Hex. impl.	Tornillo Screw	Ref.
⊕ 2,40 mm	Largo • Long	KDS OF 0007
	Corto • Short	KDS OF 0012

Tomas impresión cubeta abierta (CA) plataforma Ø4,10 mm Open tray (OT) impression coping for Ø4.10 mm platform

Hex. impl.	Tornillo Screw	Ref.
⊕ 2,70 mm	Largo • Long	KDS OF 0002
	Corto • Short	KDS OF 0009

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Tomas impresión cubeta cerrada (CC)
Closed tray (CT) impression coping

Cantidad (Ud.) Quantity (Un.)	Ref.
5	KDS OF 0014

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To replace a dental crown with a restoration dental prosthesis, a castable is used, which is defined as a meltable plastic machined attachment that, once positioned over the coronal part of the external connection implant, will be the initial basis to model the form of the final structure, to be cast in the corresponding metal.

Key features

- Machined castables in white meltable plastic material.
- Physiologically inert and excellent dimensional stability.
- Available in rotary and non-rotary version.

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Clinic screws are defined as mechanical devices used for the final anchoring of the abutments, the transfer or the castable directly to the implant, while laboratory screws are the mechanical devices used for the temporary anchoring of castables to the replica of the implant used by the dental lab where the construction of dental crowns or prostheses takes place.

Key features

- Machined grade V titanium clinic screws according to ISO-5832.
- Machined laboratory screws in stainless steel quality.
- Tightening or screws with b/f 1.20 mm hex wrench.
- Apical bit to facilitate the start of threading.
- Laboratory screw with groove to differ it from clinic screw.

Tornillos clínica Clinic screws		
Hex. llave (mm) Hex. wrench (mm)	Rosca (mm) Thread (mm)	Ref.
◇ 1,20 mm	M2,0x0,40	KDT 0C 0030
	M2,5x0,45	KDT 0C 0037
	M1,8x0,35	KDT 0C 0045

Tornillos laboratorio Laboratory screws		
Hex. llave (mm) Hex. wrench (mm)	Rosca (mm) Thread (mm)	Ref.
◇ 1,20 mm	M2,0x0,40	KDT 0C 0004
	M2,5x0,45	KDT 0C 0038
	M1,8x0,35	KDT 0C 0050

Tornillos clínica ME ME clinic screws	
Plataforma Platform	Ref.
ø4,80	KDT 0C 0056
ø6,00	KDT 0C 0060

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gmi phoenix surgical kit consists of an autoclavable case made of high strength resin, containing all the necessary material for the correct positioning of gmi phoenix external connection dental implant system. The components are arranged in such a way as to be easy to find, at all times, the necessary instrument at any time for each stage of the intervention.

gmi surgical kits are distributed in cartons, made of high strength technical polymer, and should be steam autoclaved before use at a maximum temperature of 134°C.

KYC OF 3009


	Descr.	Un.	Ref.
Fresa lanceolada	Lance-shaped drill	1	KYF OC 1224
Fresa piloto ø2,00 mm	Pilot drill ø2.00 mm	1	KYF OC 2221
Fresa guía ø2,00 / ø2,70 mm	Guide drill ø2.00 / ø2.70 mm	1	KYF OC 3227
Fresa final ø2,80 mm	Twist drill ø2.80 mm	1	KYF OC 0128
Fresa final ø3,00 mm	Twist drill ø3.00 mm	1	KYF OC 0130
Fresa final ø3,25 mm	Twist drill ø3.25 mm	1	KYF OC 0132
Fresa final ø3,50 mm	Twist drill ø3.50 mm	1	KYF OC 0135
Fresa final ø4,25 mm	Twist drill ø4.25 mm	1	KYF OC 0143
Fresa final ø5,10 mm	Twist drill ø5.10 mm	1	KYF OC 0151
Fresa final ø5,40 mm	Twist drill ø5.40 mm	1	KYF OC 0154
Llave transportador plana 45°	Manual wrench 45°	1	KYL OC 0093
Llave carraca IP	TI ratchet wrench	1	KYL OF 0113
Llave-transportador HEX - 3,00 mm IP corta	Short HEX-3.00 mm implant carrier TI wrench	1	KYL OF 0123
Llave-transportador HEX - 3,00 mm IP larga	Long HEX-3.00 mm implant carrier TI wrench	1	KYL OF 0124
Llave-transportador HEX - 3,00 mm CA	HEX-3.00 mm implant carrier HP wrench	1	KYL OF 0122
Llave fija HEX - 1,20 mm manual	HEX-1.20 mm fixed manual wrench	1	KYL OC 0011
Llave carraca corta HEX - 1,20 mm	Short HEX-1.20 mm ratchet wrench	1	KYL OF 0128
Llave carraca larga HEX - 1,20 mm	Long HEX-1.20 mm ratchet wrench	1	KYL OF 0129
Llave fija HEX - 0,80 mm manual	HEX-0.80 mm fixed manual wrench	1	KYL OC 0092
Alargador llaves CA	HP tip extender	1	KYL OC 0023
Paralelizador ø2,00 / ø2,80 mm	Paralleling pin ø2.00 / ø2.80 mm	3	KYL OC 0078
Medidor profundidad ø2,80 mm	Depth gauge ø2.80 mm	1	KYL OC 0094
Medidor profundidad ø3,50 mm	Depth gauge ø3.50 mm	1	KYL OC 0095

phoenix

implante dental de conexión externa
external connection dental implant

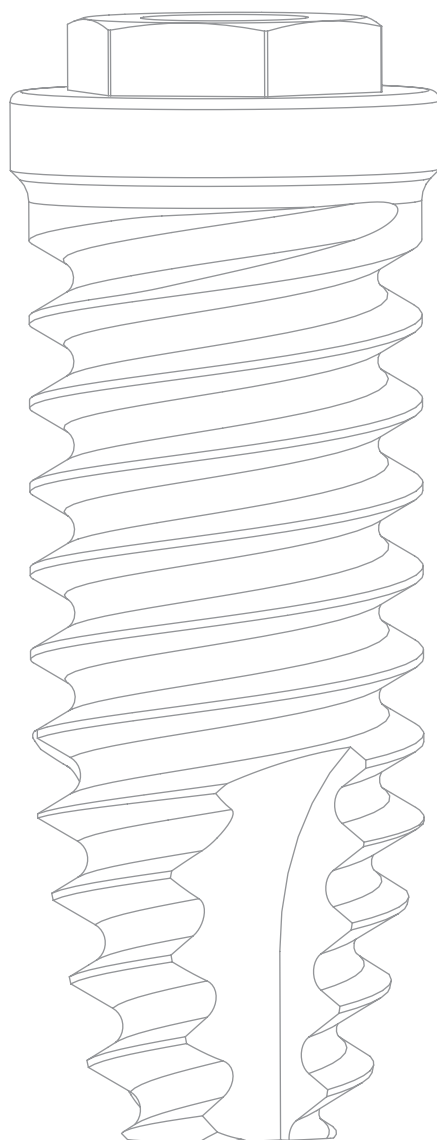Distribuido por:
Distributed by:

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Pol. Ind. El Segre
C/ Enginyer Mies 705-B
25191 Lleida (SPAIN)
Tel. (+34) 973 184 350
Fax (+34) 973 183 278
www.ilerimplant.com